

M482 & M483 HISTOSETTE® II CASSETTES IN E-Z LOAD™ SLEEVES FOR AUTOMATED LABELING INSTRUMENTS

Made of acetal polymer

This newly patented stack of cassettes allows you to load the Micro Writer® and Shur/Mark® Cassette Labeling Instruments in one simple operation. Save time and money with these convenient stacks of 50 unique **Simport** cassettes. No more manual insertion, one cassette at a time. Just load the stack in the hopper, cut and remove the holding tie and you are ready for printing.

These specially designed cassettes belong to the world-wide known **Simport** Histosette® II design. Molded from a special high density polymer, these cassettes keep specimens safely submerged in liquid and are totally resistant to the chemical action of histological solvents. The efficient flow-through slots maximize fluid exchange and ensure proper drainage.


Two types of cassettes are offered in order to suit your particular needs: a regular tissue cassette and a biopsy model designed to hold small biopsy samples securely during the embedding process.

Anterior printing area is at a 45° angle, offering an unobstructed view of the writing surface and making the cassette perfectly suitable to be used with cassette labeling instruments.

Two types of lids are available. The one-piece disposable plastic cover eliminates the need for reusable steel lids. It can be opened and closed as often as necessary and it always relocks securely without danger of specimen loss. But for those who prefer a reusable metal lid, please refer to M481 model on page 17. Available in non-cytotoxic, non-metallic popular colors.

Each case contains two boxes of 20 stacks.

Order cassette bases and covers separately.

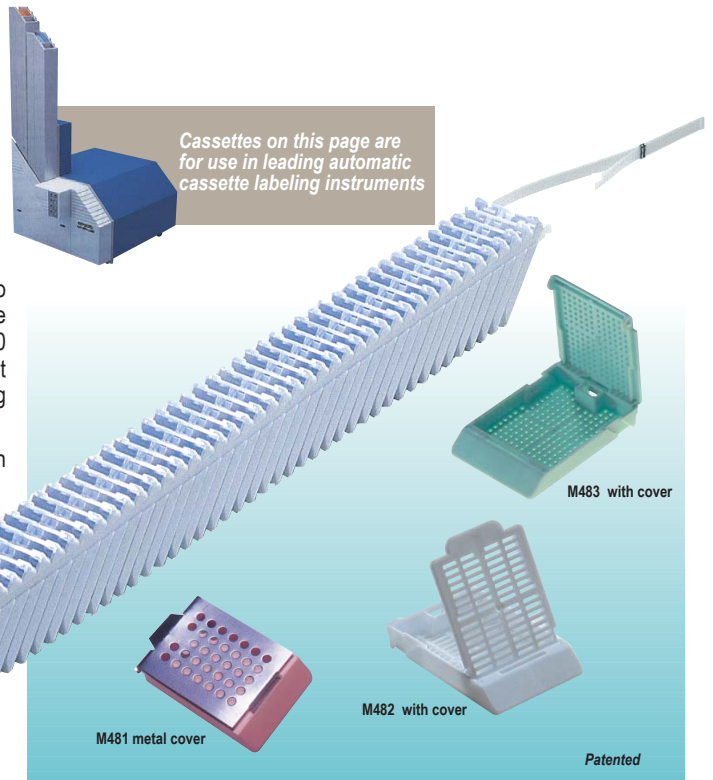

TISSUE CASSETTES (In stacks of 50)				BIOPSY CASSETTES (In stacks of 50)			
Base Cat. #	Cover Cat. #	Color	Qty/Cs	Base Cat. #	Cover Cat. #	Color	Qty/Cs
M482-2BA	M482-2LI	White	2000	M483-2BA	M483-2LI	White	2000
M482-3BA	M482-3LI	Pink	2000	M483-3BA	M483-3LI	Pink	2000
M482-4BA	M482-4LI	Green	2000	M483-4BA	M483-4LI	Green	2000
M482-5BA	M482-5LI	Yellow	2000	M483-5BA	M483-5LI	Yellow	2000
M482-6BA	M482-6LI	Blue	2000	M483-6BA	M483-6LI	Blue	2000
M482-7BA	M482-7LI	Peach	2000	M483-7BA	M483-7LI	Peach	2000
M482-8BA	M482-8LI	Tan	2000	M483-8BA	M483-8LI	Tan	2000
M482-9BA	M482-9LI	Gray	2000	M483-9BA	M483-9LI	Gray	2000
M482-10BA	M482-10LI	Lilac	2000	M483-10BA	M483-10LI	Lilac	2000
M482-11BA	M482-11LI	Orange	2000	M483-11BA	M483-11LI	Orange	2000
M482-12BA	M482-12LI	Aqua	2000	M483-12BA	M483-12LI	Aqua	2000

How to insert HISTOSETTE® II cassettes in EZ-LOAD™ sleeves

1 IMPORTANT

For magazines with metal release mechanism
LOAD CASSETTE SLEEVE FROM BOTTOM OF MAGAZINE

Bottom of magazine

2

Cut plastic tie, lift and discard while holding cassettes in place with forefinger.

3

Place hopper in labeling instrument.

How to handle HISTOSETTE® II cassettes after printing

HOW TO INSTALL COVER

Insert vertically until hinges click firmly into slots in front of cassette.

Cover can also be inserted at an angle.

Press down on back tab.

HOW TO REMOVE COVER

Lift tab at the back of cassette with forefinger while depressing the center of lid with thumb.

To remove the cover, pivot the lid forward and it will disconnect automatically.